

## Facial Steamer Operating Instructions



Thank you for purchasing NTFS-618 facial steamer , please read carefully the instructions before you use ,to fully understand the method of use.

Hydrating and skin surface facial steamer NTFS-618  
Product material: environmental protection plastic ABS / PP  
/high performance PTC  
Product specifications : maximum 50ml

- 1.For pore cleansing,deep hydration&purification of the face
- 2.0/1/2 Switch setting with LED indicator light
- 3.Wide facial mask uses gentle steam to open clogged pores and rejuvenate skin
- 4.Fashion appearance & easy to use
- 5.Overheating protection device, prevent the dry burn

### Precautions

- 1.If you suffer from heart disease, severe rosacea, eczema or psoriasis patients do not use this product
- 2.Pregnant women and patients with serious diseases, such as in the use of the products need to add herbs and essential oil before, please consult a physician.
- 3.Once appear the skin inflammation or discomfort, immediately stop using essential oils and herbs.
- 4.This product only suitable for family use, not in wet environment use.
- 5.The product before turning on the power supply voltage and check the family product plug shows voltage is consistent
- 6.Such as power cord damage, to avoid danger must by manufacturer or its maintenance department or similar full-time personnel to change

Hydrating and skin surface facial steamer NTFS-618

Product material: environmental protection plastic ABS / PP  
/high performance PTC

Product specifications : maximum 50ml


- 1.For pore cleansing,deep hydration&purification of the face
- 2.0/1/2 Switch setting with LED indicator light
- 3.Wide facial mask uses gentle steam to open clogged pores and rejuvenate skin
- 4.Fashion appearance & easy to use
- 5.Overheating protection device, prevent the dry burn

## **Precautions**

- 1.If you suffer from heart disease, severe rosacea, eczema or psoriasis patients do not use this product
- 2.Pregnant women and patients with serious diseases, such as in the use of the products need to add herbs and essential oil before, please consult a physician.
- 3.Once appear the skin inflammation or discomfort, immediately stop using essential oils and herbs.
- 4.This product only suitable for family use, not in wet environment use.
- 5.The product before turning on the power supply voltage and check the family product plug shows voltage is consistent
- 6.Such as power cord damage, to avoid danger must by manufacturer or its maintenance department or similar full-time personnel to change

7. Non-professionals do not open this product, avoid life safety.
8. When the steam has completely evaporated or once you are finished, turn the switch to off, then unplug from the electrical outlet

## Component

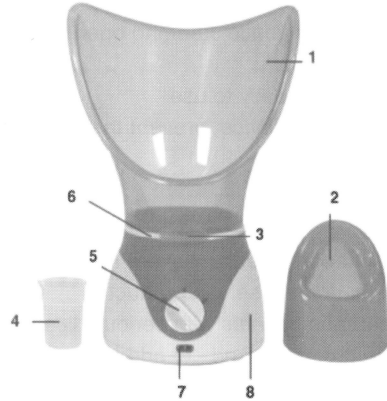

- |                  |                          |
|------------------|--------------------------|
| 1. Face Mask     | 5. Switch intorum        |
| 2. Nasal Mask    | 6. Stand                 |
| 3. Valve         | 7. Power indicator light |
| 4. Measuring cup | 8. Base                  |

## Product Efficacy

1. For pore cleansing, deep hydration & purification of the face
2. Dredge pores, remove grease in pores, deep clean the skin
3. Promote the face blood circulation, make skin renew

4. Replenishing water, maintain water oil balance of skin, increase skin elasticity
5. Relieve a cold, sinus problems
6. Can use with incense or essential oils help to eliminate fatigue, eliminate pouch and black rim of the eye, in the home can enjoy aromatherapy sauna

## Method of use

1. Load 50 ml of water into the base, Do not fill water above the MAX line show inside the fill area of the base
2. The bottom of the steam mask cartoon aligned at the bottom of the slot, rotate clockwise fixed.
3. When using a nasal mask, the mask can be placed directly on the middle, you can adjust the steam direction
4. Turning on the power, you can choose 1 or 2 switch setting, 1 low, 2 high
5. When the steam has completely evaporated or once you are finished, turn the switch to off, then unplug from the electrical outlet
6. Relax for 8 minutes after the hydration process, and then clean your face with cold water

## Product performance

1. Product specifications: 110-240V 50/60Hz 130w
2. PTC heating
3. CE / ROHS certification
4. 0/1/2 Switch to adjust the amount of steam-conditioning
5. Overheating protection device, prevent the dry burn
6. Adjustable steam size, steam choice more free and flexible

BE SURE TO READ ALL INSTRUCTIONS BEFORE USING.  
WARNING: DO NOT LEAVE UNATTENDED, KEEP FROM CHILDREN.